



# VISIRAY

## Radiopaque & Phosphorescent Rulers

Radiopaque and phosphorescent rulers for quality control of congruence between light, radiation and x-ray fields.

### Features:

- Radiopaque markings (1mm)
- Visible markings (1mm)
- Phosphor coating (28-120 kV)
- Removable protective casing

### Specification:

- Size: 45cm, 30cm, 20cm and 10cm
- Bespoke sizes on request

### Note:

- The phosphor exhibits no significant afterglow
- The room must be darkened
- A video capture device is required (either our 4K camera, or your own camera e.g. a mobile phone)





# VISIRAY

## Radiopaque & Phosphorescent Rulers




Protective cases can be removed by removing the magnetically attached cover and extracting the ruler.



# VISIRAY

Radiopaque & Phosphorescent Rulers

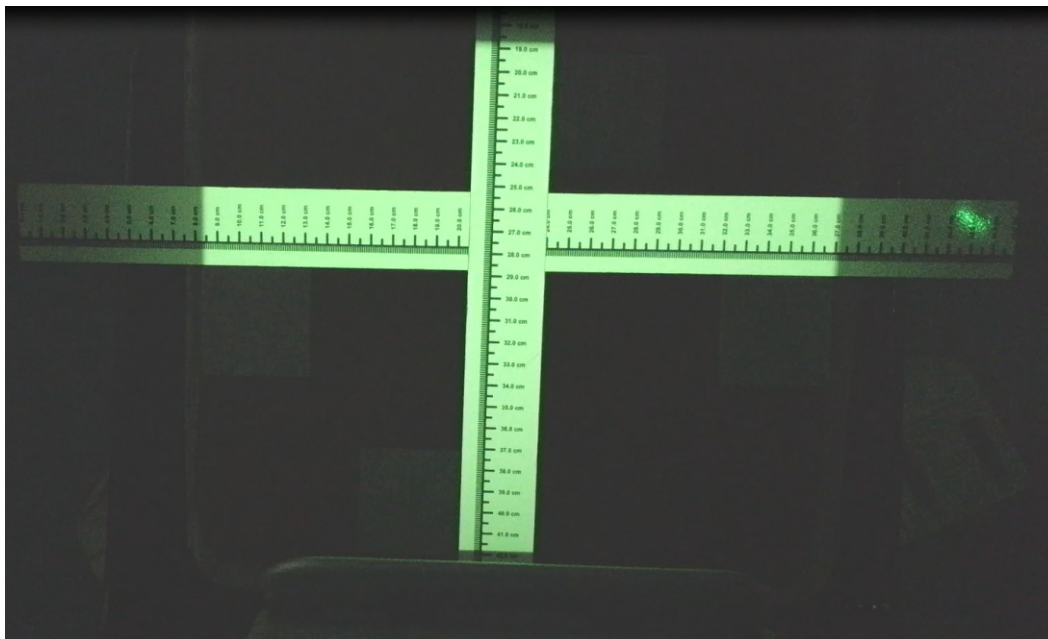

Frame captured from 4K resolution camera  
showing radiation field size




**4K Camera  
Available**



# VISIRAY

## Radiopaque & Phosphorescent Rulers

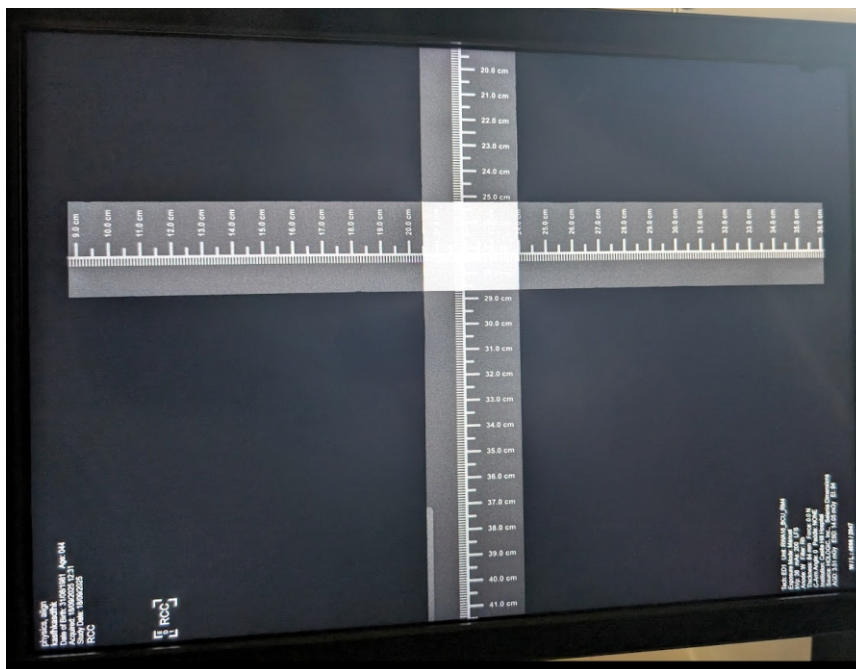

Radiograph showing x-ray image field size